

DEBUBBLERS & DEGASSERS

**ACTIVE
DEBUBBLERS**
PAGE 175

**APPLICATION
NOTE — WHY
DEGAS YOUR
MOBILE PHASE?**
PAGE 177

**STAND ALONE
DEGASSERS**
PAGE 178



How to Order

Direct from IDEX Health & Science

For North America orders:

Tel: +1 800 426 0191 / +1 360 679 2528

Fax: +1 800 359 3460 / +1 360 679 3830

E-mail: CustomerService.hs@idexcorp.com

For Europe orders:

Tel: +49 (0) 1801 808 800

Fax: +49 (0) 9377 1388

E-mail: CustomerService.hsEurope@idexcorp.com

For Asia orders:

E-mail: CustomerService.hsAsia@idexcorp.com

NORTH ASIA

Tel: +86 10 6566 9090 / +86 21 5241 5599

Fax: +86 10 6567 5152 / +86 21 5241 8338

SOUTHEAST ASIA AND OCEANIA

Tel: +65 6684 7307

Fax: +65 6795 8070

SOUTH ASIA

Tel: +91 22 66 78 55 50

Fax: +91 22 66 78 00 55

JAPAN

Tel: +81 48 240 5750

Fax: +81 48 259 0715

For Online orders:

Web: www.idex-hs.com

Please include the following information with your order: Purchase Order Number, Contact Name, Shipping and Billing Address, Telephone Number, Product Number(s) and Quantity(ies). If you are a Purchasing Agent, please include the name of the person for whom the products are ordered.

For your information, our U.S. Federal Taxpayer Identification number is 01-0736657; European Taxpayer Identification number is DE 146581180.

Order Processing/Shipping

Orders process within one business day and ship within two business days. Requests for earlier shipments can be made — please contact IDEX Health & Science for more details.

All shipments are FCA unless otherwise noted. Shipping charges will be prepaid and added to the invoice. Or, if requested, we will charge your carrier account number.

Terms

Terms of payment are net 30 days. **Credit cards are accepted for orders placed directly with our U.S. site. No international credit card payment is possible.** Customers are responsible for C.O.D. charges if they elect this payment method. Terms and prices are subject to change without notice. Visit www.idex-hs.com for complete Terms & Conditions.

Return of Materials

Returns must be authorized in advance. Please contact us within 30 days of purchase for your Returns Material Authorization (RMA) Number. A restocking fee may apply.



Customer Number:

Direct from our Global Distribution Network

The IDEX Health & Science brands represented in this catalog (Ismatec®, Isolation Technologies™, Rheodyne®, Systec®, and Upchurch Scientific®) have established relationships with a large, global network of Distributors, many of whom can offer the product support and assistance that we cannot, including:

- ▶ Carrying local inventory of commonly-ordered items
- ▶ Understanding YOUR working environment and requirements
- ▶ Providing technical application assistance
- ▶ Answering questions regarding new and existing products
- ▶ Complete availability during your time zone business hours
- ▶ Personal visits as required
- ▶ Access to complementary product lines to provide a "one-stop shopping" experience

We have come to rely on our authorized Distributors to provide value-added service to end-users of our products. So, while we are happy to take your order directly, we encourage you to contact your local Distributor — we feel confident it will be a positive experience!

How to Find A Local Distributor

For a complete distributor listing please visit www.idex-hs.com/distributors or email:

US: CustomerService.hs@idexcorp.com

Europe: CustomerService.hsEurope@idexcorp.com

Asia: CustomerService.hsAsia@idexcorp.com

Our customer service team will help you contact a distributor in your area.

Return of Materials

Product purchased through a Distributor must be returned directly to that Distributor. Please contact them for Return Material Authorization procedures.

The products in this catalog are intended for use with analytical, biotechnology and diagnostic equipment and accessories. THE MANUFACTURERS OF THE PRODUCTS IN THIS CATALOG SHALL HAVE NO LIABILITY WHATSOEVER DUE TO ANY MISUSE OF THEIR PRODUCTS. Of course, the safe use of our products depends on our customers, since it is you who select and control the protective gear and safety procedures used, as well as the pressures, temperatures, solvents, samples, ventilation, and other variables. Product and material performance ratings are provided as guides only. Individual field tests should be performed by customers to determine safe operating parameters given your particular procedures and use.

IDEX Health & Science facilities are certified ISO 9001. Some facilities are also certified ISO 13485:2003.



Visit us online!
www.idex-hs.com

Systec® Debubbler Series

Remove Bubbles, Dissolved Gas, or Both!

Dissolved gases and bubbles in system liquids cause dispense volume anomalies in many instruments, negatively affecting both dispense precision and analytical accuracy. Now you have a choice of components for actively removing bubbles with or without also removing dissolved system gases. Online Vacuum Degassing offers operating convenience, high efficiency and low operating costs compared to other common degassing technologies.

Debubbler/Degasser

Combines Vacuum Degassing with Active Bubble Removal

- ▶ Improves instrument performance — reduces downtime due to bubble formation.
- ▶ Fewer false positives due to reduction of partial reagent dispenses.
- ▶ Easily integrates into any pump, degassing tray, or stand-alone degassing application.

Active Debubbler


Remove Bubbles in Fluid Stream Before or After the Pump

- ▶ Improves instrument performance — reduces downtime due to bubble formation.
- ▶ Fewer false positives due to reduction of partial reagent dispenses.
- ▶ Easily integrates into any pump, degassing tray, or stand-alone degassing application.

Transfer-Line Degasser

Removes Dissolved Gases During Fluid Transfer

- ▶ Eliminates baseline fluctuations for improved detector sensitivity.
- ▶ Coaxial design reduces number of connections, improves reliability.
- ▶ Single lumen design increases degassing reliability.

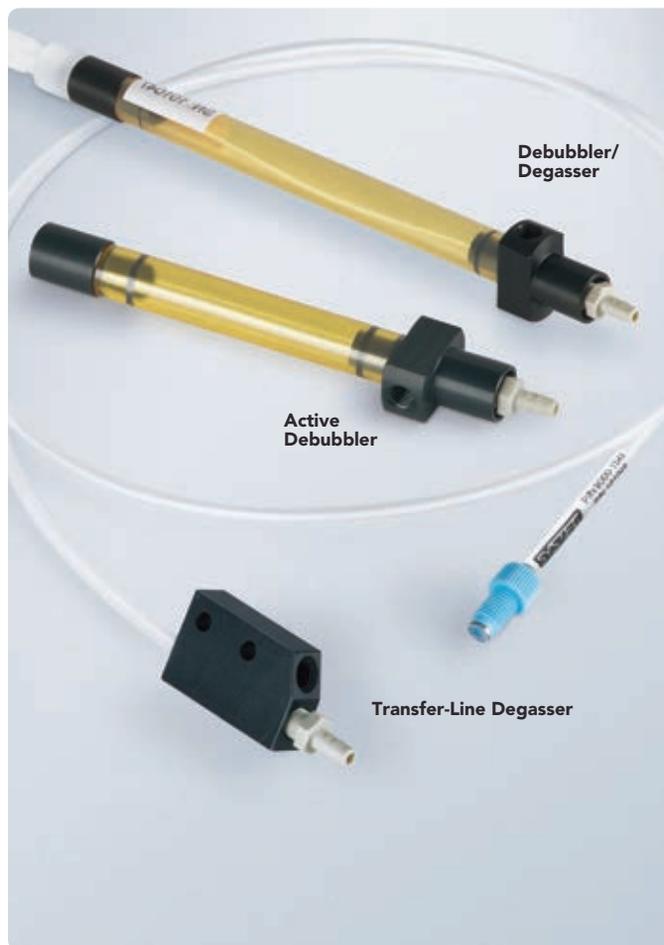


APPLICATION NOTE

- ▶ Liquid handling
- ▶ IVD
- ▶ HPLC/UHPLC
- ▶ O₂ and CO₂ removal

In medical analyzers, bubbles interfere with critical volumetric reagent dispenses and cause sample failures, wasting time and money. Because bubbles adhere to nearly every part of a dispensing system, high velocity or induced turbulent flow is often used to displace and discharge bubbles from the flow stream and into a waste area. These alternative processes waste reagents and are time consuming, unpredictable, and may additionally require designing the system to recognize bubbles are present. Regardless of how the systems are designed, aqueous systems will always be subject to the laws of physics that cause out-gassing during changes in fluid temperature, pressure, or chemicals mixture. In fluid applications like these, debubblers are the optimal solution to capture and remove formed bubbles to prevent sample dispense inaccuracies, and degassing is ideal to prevent downstream bubble formation from recurring.


Part No.	Description	Standard Bubble Trap Size	Transfer Line Length	Internal Volume	Max Bubble Capacity
DEBUBLER SERIES – AVAILABLE STANDARD CONFIGURATION					
9000-1540	2.5 mL Active Debubbler	2.5 mL	—	2.5 mL	2.5 mL
9000-1541	5 mL Active Debubbler	5 mL	—	5 mL	5 mL
9000-1544	2.5 mL Debubbler/Degasser	2.5 mL	17.5" (444.5 mm)	2.5 mL in transfer line + 2.5 mL in bubble trap	2.5 mL
9000-1545	5 mL Debubbler/Degasser	5 mL	34" (863.6 mm)	5 mL in transfer line + 5 mL in bubble trap	5 mL
★ 9000-1549	1.1 m Transfer-Line Degasser	—	1.1 m (43")	4 mL	N/A

Overall Dimensions

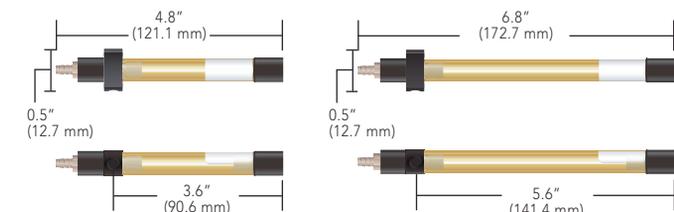


9000-1540
Active Debubbler
2.5 mL


9000-1541
Active Debubbler
5 mL

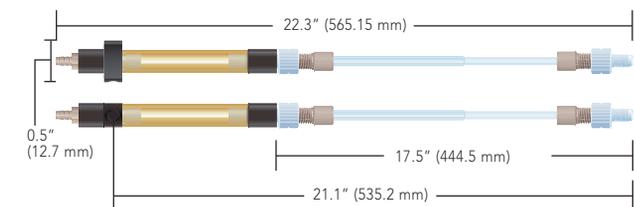

9000-1544
Debubbler / Degasser
2.5 mL

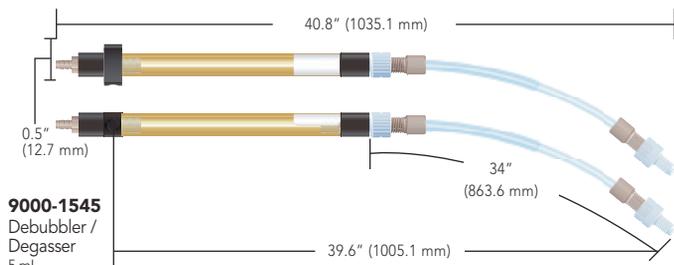

9000-1545
Debubbler / Degasser
5 mL

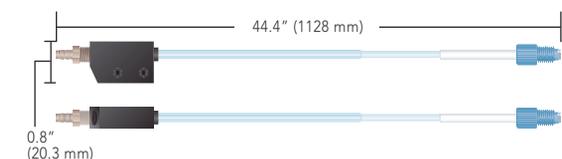

9000-1549
Transfer-Line Debubbler
1.1 meter


Please note: These drawings are not actual size.

BENEFITS


	Active Debubbler	Debubbler/ Degasser	Transfer-Line Degasser
Perfect for applications that require dissolved gas like oxygen for reaction kinetics	✓		
Improves dispense precision by capturing and removing bubbles	✓	✓	
Eliminates false positives and reduces reagent waste by improving instrument performance	✓	✓	
Easily integrates into fluidic path	✓	✓	✓
Creates stable instrument performance across system and environmental fluctuations	✓	✓	✓
Prevents the formation of bubbles downstream of the degasser		✓	✓
Eliminates fluctuations for improved detector sensitivity and accuracy by preventing bubble formation		✓	✓
Improves fluidic system reliability because coaxial design reduces the number of connections			✓
Flexible design can be implemented as transfer line in new instruments or existing instruments that don't have space available			✓
Minimizes fluidic system internal volumes to reduce reagent cost			✓

IMPLEMENTATIONS


Typical Debubbler Implementation

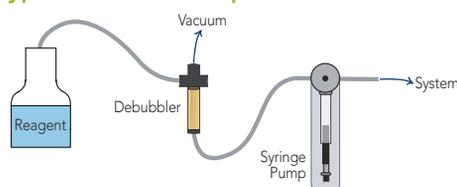

Transfer-Line Degasser Implementation

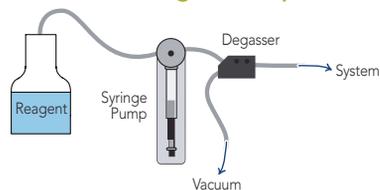

Gas bubbles are actively removed from a flowing liquid stream by vacuum via the PORIDEX membrane.

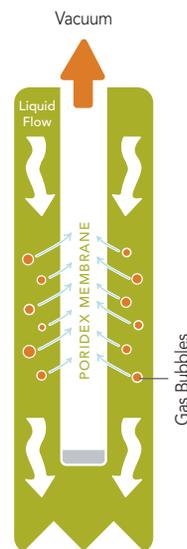

NOTE

Degassing tubing is flexible and therefore can be coiled to shorten the overall length or used to transfer the fluid within an instrument to the next desired location.

SPECIFICATIONS (ALL PLATFORMS)

	Active Debubblers	Degasser/ Debubblers	Transfer-Line Degasser
Bubble Removal (volume of air removed/min @ 10 mL/min H₂O)	Up to 30 cc	Up to 30 cc	N/A
Degassing Efficiency† @ 1 mL/min H₂O	N/A	2.5 mL size: 36% O ₂ removal 5.0 mL size: 55% O ₂ removal	< 4 ppm dissolved O ₂ at 5 mL/min
Membrane Material	PORIDEX®	PORIDEX	PORIDEX
Wetted Materials	PORIDEX, Polyolefin, FEP, ETFE, Ultem®	PORIDEX, Polyolefin, FEP, ETFE, Ultem	PORIDEX, Polyolefin, FEP, ETFE
Solvent Compatibility	Solutions > 50% aqueous. Not compatible with detergent concentrations > 0.05%.		
Standard Bubble Trap Volume	2.5 / 5.0 mL	2.5 / 5.0 mL	N/A
Transfer-Line Volume	N/A	2.5 / 5.0 mL	< 4 mL
Maximum Operating Pressure	200 kPa (30 psi) @ 25 °C		
Maximum Operating Temperature	40 °C		
Recommended Vacuum Level	Minimum 16 kPa absolute		
Liquid Connection	1/4-28 fitting system		
Vacuum Connection	Tubing vacuum port(s) for 1/8" (3 mm) ID elastomeric tubing		
Pressure Drop	0.8 mm Hg / mL / min (assumes laminar flow and viscosity of 1 cP)		

† Degubbling / degassing efficiency can be optimized based on flow rate, fluid to be degassed, and gas to be removed.

APPLICATION NOTE

Why Degas Your Mobile Phase?

Dissolved air in HPLC mobile phases can result in flow rate instability and baseline disturbance.

Flow rate instability: Non-degassed mobile phase can outgas in the pump head, causing bubbles to be formed and trapped inside the head or check valves. These bubbles can cause flow disturbances and pressure fluctuations, resulting in flow rate instability.

Baseline disturbance: As the mobile phase passes through the column, it experiences a large pressure drop. Non-degassed mobile phase can outgas due to this pressure differential, causing air bubbles to form. Air bubbles passing through or lodging in the flow cell cause detection disturbances, exhibited as baseline noise.

Why Use a Degassing System?

Helium sparging is a common means of degassing HPLC solvents. This method has its drawbacks, however. Helium tanks are expensive and bulky, and solvent backup and contamination are concerns. In addition, helium sparging can change the composition of a premixed mobile phase over time, due to the difference in the evaporation rates of mobile phase components.

In contrast, the Systec® Degassing System has none of these drawbacks, and it is extremely fast and efficient at removing dissolved gases — more efficient than helium sparging or PTFE-based degassing systems.



Tubing Connections

We recommend ETFE tubing (page 73) be used to limit regassing of mobile phase between the degasser and your pump. ETFE is recommended because of its superior impermeability to gases (compared to PTFE, FEP, and PFA tubing). Applicable flangeless fittings for 1/8" OD tubing are found on page 26.

GPC and HFIP Applications

Standard degassing chambers, with PEEK bulkhead unions, are not recommended for GPC applications or for use with HFIP (hexafluoroisopropanol). Special GPC "hardened" versions are available. Please contact us for more information.

TYPICAL DEGASSER IMPLEMENTATION

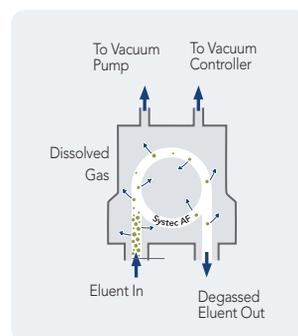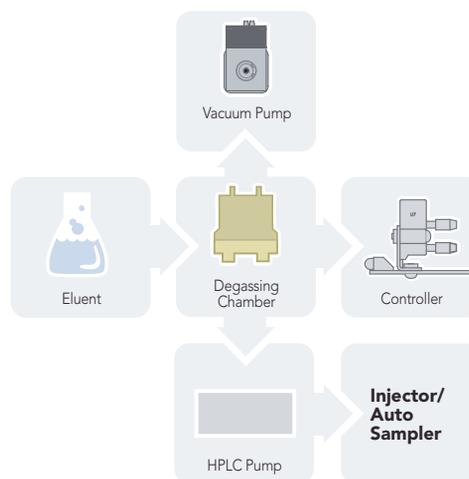

Dissolved gases are actively removed from a flowing liquid stream by vacuum via the Systec AF® membrane.

System[®] Stand Alone MINI & Prep Scale Vacuum Degassing Systems

- ▶ Analytical and Prep scale models
- ▶ Ultra-high degassing efficiency
- ▶ Low volume, easy to prime
- ▶ Patented control eliminates baseline fluctuations
- ▶ Inert flow path
- ▶ 5+ year lifetime

The System Stand-Alone MINI and Prep-Scale HPLC vacuum degassing systems are high-efficiency, in-line modules that remove dissolved gases from the mobile phase. Their unique design assures reliable continuous operation and the highest level of performance available without the need for helium sparging. Up to five solvent lines may be degassed simultaneously by one unit.

ZHCR[®] Control with Built-in Test Diagnostics

- ▶ Microcontroller self-test vacuum sensor validation on power-up
- ▶ Continuous vacuum system monitoring to ensure optimum operational conditions are maintained
- ▶ Vacuum system fault detection and shutdown function indicators

AF / ZHCR Degassing Technology


Flow-through vacuum degassing chamber with a single amorphous perfluorinated copolymer (System AF[®]) degassing membrane, enabling degassing efficiency 50 times that of PTFE.

The ZHCR (Zero Hysteresis / Constant Run) vacuum pump employs a patented closed-loop, micro-stepping rpm control strategy permitting the pump to run with continuously variable speed, providing quick pull-down at high rpm, and then sustaining a consistent vacuum level at low rpm.

Fluctuations in detector baseline due to changes in vacuum level are eliminated by not having to repeatedly stop and start a single-speed pump. This also greatly reduces wear and noise.

The brushless motor enables quiet operation and is appropriate for environments where solvent vapors may be present.

Overall Dimensions

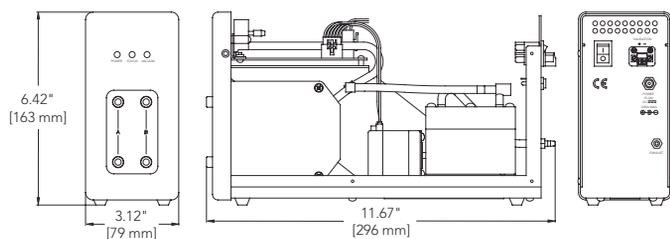

SPECIFICATIONS

Stand-Alone Degassing Modules	
Maximum Number of Degassing Channels	5
Degassing Efficiency [†] @ 1 mL / min MeOH	> 70% O ₂ Removal
Membrane Material	SYSTEC AF
Other Wetted Materials	PEEK, PPS(GF), PTFE(GF), FEP
Solvent Compatibility	Not compatible with fluorinated solvents. Special version available for GPC solvents.
Flow Path ID	1.14 mm (0.045")
Internal Volume	480 µL (standard)
Maximum Pressure (@ 25 °C)	0.5 MPa (70 psi)
Pressure Drop	0.18 kPa/mL/min

[†] Degassing efficiency can be optimized based on flow rate, fluid to be degassed, and gas to be removed.
¹ Standard ID; other sizes available.


NOTE

Power Requirement

Input Power required with AC Adapter (included): 100 to 240 V AC (±10%), 1A, 50 to 60 Hz (±3 Hz). Four interchangeable wall sockets are supplied with the AC Adapter: North America/Japan, U.K., Continental Europe, Australia.

CE Certification

This product has been certified under the following CE testing standards: EN61326-1; EN55011; EN61300-3-2; EN61300-3-3, & EN61010-1.

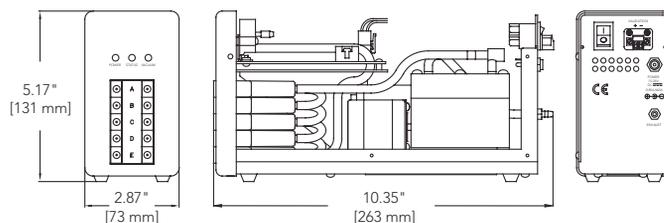

Part No.	Number of Channels	Channel Volume	Max HPLC Gradient Flow Capability	Pressure Drop ^E	Degassing Flow Path ID
SYSTEC STAND ALONE MINI VACUUM DEGASSING SYSTEMS — AVAILABLE CONFIGURATIONS^A					
★ 0001-6500	2	480 µL	2.0 mL/min ^C	0.18 kPa/mL/min	0.045" (1.14 mm)
★ 0001-6501	4	480 µL	2.0 mL/min ^C	0.18 kPa/mL/min	0.045" (1.14 mm)
SYSTEC STAND ALONE PREP SCALE VACUUM DEGASSING SYSTEMS — AVAILABLE CONFIGURATIONS^{A, B}					
0001-6482	2	8.4 mL	20 mL/min ^D	0.04 kPa/mL/min	0.065" (1.65 mm)
★ 0001-6484	2	13.8 mL	40 mL/min ^D	0.06 kPa/mL/min	0.065" (1.65 mm)

A. Custom configurations are available. Consult us for your own OEM solution to your specific application.

B. The standard prep scale chambers are not recommended for GPC applications or for use with HFIP (hexafluoroisopropanol).

C. The flow rates given are for a gradient mixture of 50/50 MeOH/H₂O, with a typical low pressure gradient mixing valve. Higher flow rates are possible with high pressure mixing.

D. The flow rates given are for a gradient mixture of 60/40 MeOH/H₂O, with a typical low pressure gradient mixing valve. Higher flow rates are possible with high pressure mixing.

E. Estimated tubing pressure per unit change in flow assuming laminar flow with a viscosity of 1.0 cP